


TITRANT H+ /M = 0.05019

| V/ML | E/MV | PH | NBAR | TH | TH-NBAR |
|------|--------|------|------|------|---------|
| 0.20 | -220.4 | 9.34 | 0.07 | 0.09 | 0.01 |
| 0.40 | -205.3 | 9.08 | 0.14 | 0.15 | 0.01 |
| 0.60 | -194.9 | 8.90 | 0.21 | 0.22 | 0.01 |
| 0.80 | -186.8 | 8.77 | 0.28 | 0.29 | 0.01 |
| 1.00 | -180.3 | 8.66 | 0.35 | 0.35 | 0.01 |
| 1.20 | -174.9 | 8.57 | 0.41 | 0.42 | 0.00 |
| 1.40 | -169.2 | 8.47 | 0.48 | 0.49 | 0.01 |
| 1.60 | -164.4 | 8.39 | 0.55 | 0.56 | 0.01 |
| 1.80 | -160.0 | 8.31 | 0.62 | 0.63 | 0.01 |
| 2.00 | -155.7 | 8.24 | 0.69 | 0.70 | 0.01 |
| 2.20 | -151.9 | 8.18 | 0.76 | 0.77 | 0.01 |
| 2.40 | -148.1 | 8.11 | 0.83 | 0.83 | 0.01 |
| 2.60 | -144.4 | 8.05 | 0.90 | 0.90 | 0.00 |
| 2.80 | -140.9 | 7.99 | 0.97 | 0.97 | 0.00 |
| 3.00 | -136.8 | 7.92 | 1.04 | 1.04 | 0.01 |
| 3.20 | -131.1 | 7.83 | 1.10 | 1.15 | 0.04 |
| 3.40 | -129.6 | 7.80 | 1.17 | 1.17 | 0.00 |
| 3.60 | -125.0 | 7.72 | 1.24 | 1.25 | 0.01 |
| 3.80 | -121.1 | 7.66 | 1.31 | 1.32 | 0.01 |
| 4.00 | -117.0 | 7.59 | 1.38 | 1.39 | 0.01 |
| 4.20 | -112.3 | 7.51 | 1.45 | 1.46 | 0.01 |
| 4.40 | -107.3 | 7.42 | 1.52 | 1.53 | 0.01 |
| 4.60 | -101.5 | 7.33 | 1.59 | 1.60 | 0.02 |
| 4.80 | -95.0 | 7.22 | 1.66 | 1.68 | 0.02 |
| 5.00 | -86.7 | 7.08 | 1.72 | 1.75 | 0.03 |
| 5.20 | -77.3 | 6.92 | 1.79 | 1.82 | 0.03 |
| 5.40 | -65.5 | 6.72 | 1.86 | 1.88 | 0.02 |
| 5.60 | -39.2 | 6.27 | 1.93 | 1.96 | 0.02 |
| 5.80 | 87.6 | 4.13 | 1.99 | 2.00 | 0.01 |
| 6.00 | 131.6 | 3.39 | 2.00 | 2.00 | 0.00 |
| 6.20 | 147.5 | 3.12 | 2.00 | 2.00 | 0.00 |
| 6.40 | 157.9 | 2.94 | 2.00 | 2.00 | 0.00 |
| 6.60 | 164.6 | 2.85 | 2.00 | 2.00 | 0.00 |

TITRATION 2

E0 /MV = 333.3
 INIT VOL /ML = 20.00
 LIG IN CELL/MMOL = 0.09860
 TITRANT H+ /M = 0.05019

| V/ML | E/MV | PH | NBAR | TH | TH-NBAR |
|------|--------|------|------|------|---------|
| 0.20 | -211.2 | 9.20 | 0.11 | 0.12 | 0.01 |
| 0.40 | -195.4 | 8.94 | 0.21 | 0.21 | 0.00 |
| 0.60 | -183.0 | 8.73 | 0.31 | 0.31 | 0.00 |
| 0.80 | -174.3 | 8.59 | 0.41 | 0.40 | -0.01 |
| 1.00 | -166.7 | 8.45 | 0.51 | 0.51 | 0.00 |
| 1.20 | -160.8 | 8.35 | 0.61 | 0.59 | -0.02 |
| 1.40 | -154.0 | 8.24 | 0.71 | 0.71 | -0.01 |
| 1.60 | -148.3 | 8.14 | 0.81 | 0.80 | -0.01 |
| 1.80 | -142.7 | 8.05 | 0.92 | 0.91 | -0.01 |
| 2.00 | -137.7 | 7.96 | 1.02 | 1.00 | -0.02 |
| 2.20 | -132.7 | 7.88 | 1.12 | 1.09 | -0.03 |

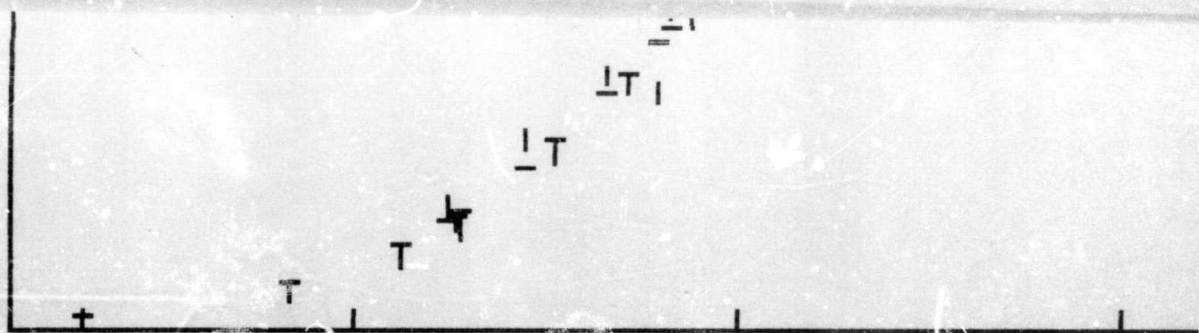

X - AXIS STARTS AT -9.9 UNITS = 1.0
 Y - AXIS STARTS AT 0 UNITS = 0.5

T = EXPERIMENTAL NBARS
 T = THEORETICAL NBARS

332.5
 19.00
 0.8733
 0.8664
 0.8591


| V/ML | E/MV | LG | NBAR | TH | TH-NBAR |
|------|------|-------|------|------|---------|
| 1.10 | 78.9 | -9.70 | 0.01 | 0.01 | 0.01 |
| 1.20 | 62.4 | -9.16 | 0.04 | 0.05 | 0.01 |
| 1.30 | 53.4 | -8.87 | 0.98 | 0.09 | 0.02 |
| 1.40 | 48.8 | -8.73 | 0.12 | 0.12 | 0.01 |
| 1.60 | 40.3 | -8.47 | 0.19 | 0.20 | 0.01 |
| 1.80 | 34.0 | -8.29 | 0.27 | 0.28 | 0.01 |
| 2.00 | 28.0 | -8.12 | 0.35 | 0.37 | 0.02 |
| 2.40 | 18.0 | -7.85 | 0.50 | 0.51 | 0.01 |
| 2.60 | 12.8 | -7.72 | 0.58 | 0.59 | 0.01 |
| 2.80 | 7.7 | -7.59 | 0.66 | 0.66 | 0.00 |
| 3.00 | 2.4 | -7.46 | 0.74 | 0.72 | -0.01 |
| 3.20 | -4.2 | -7.29 | 0.81 | 0.80 | -0.02 |

TITRATION 3

E_0 /MV = 332.1
 INIT VOL. /ML = 20.00
 H^+ IN CELL/MMOL = 0.07530
 LIG IN CELL/MMOL = 0.02910
 M IN CELL. /MMOL = 0.01650
 TITRANT B- /M = 0.05150

| V/ML | E/MV | LG | NBAR | TH | TH-NBAR |
|------|-------|-------|------|------|---------|
| 0.40 | 33.8 | -8.71 | 0.11 | 0.13 | 9.02 |
| 0.50 | 17.3 | -8.20 | 0.26 | 0.32 | 0.06 |
| 0.60 | 6.9 | -7.90 | 0.42 | 0.49 | 0.07 |
| 0.72 | -3.5 | -7.62 | 0.60 | 0.65 | 0.05 |
| 0.80 | -11.4 | -7.40 | 0.72 | 0.75 | 0.03 |
| 0.90 | -27.9 | -6.92 | 0.87 | 0.90 | 0.03 |

Log $K_1(CuII)$ -DHP-15-ane-N203

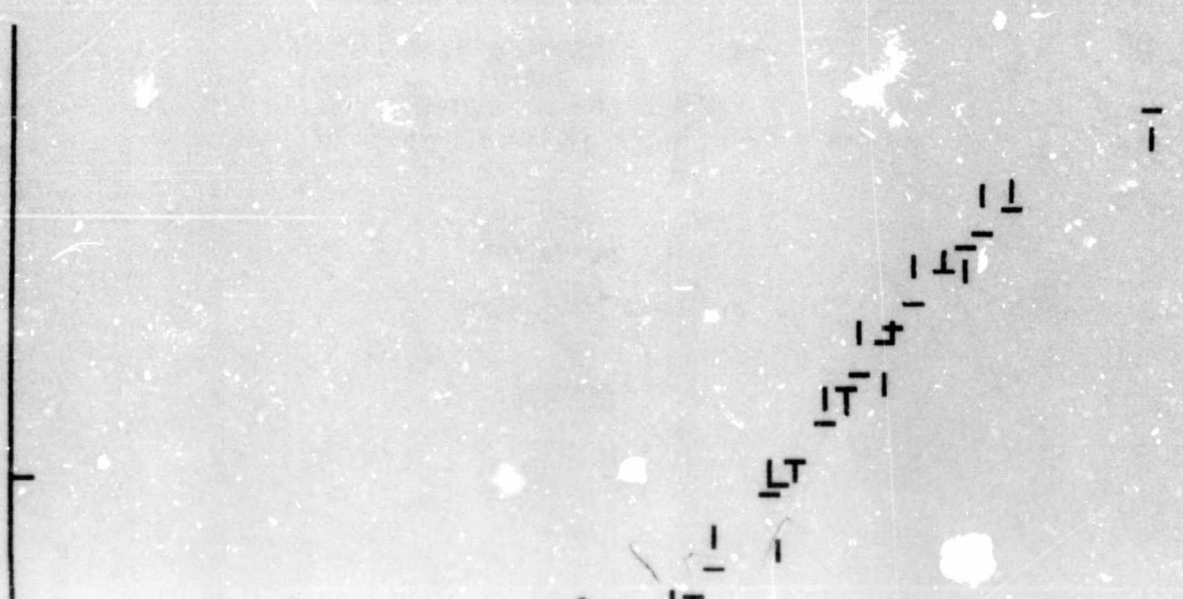

Log K1(CuII)-DHP-15-ane-N203

ACID CONSTANTS:

PKA (1) = 8.39
 PKA (2) = 7.62

MATRIX CONDITION NUMBER = 1.0000
 CONVERGENCE IN *3* CYCLES
 RMSD = .0305
 SUM DELTA SQUARED = .03

PARAMETERS AND ERRORS

LOG K (1) = 7.879 +- .013

TITRATION 1

E0 /MV = 332.5
 INIT VOL /ML = 19.00
 LIG IN CELL/MMOL = 0.07395
 M IN CELL /MMOL = 0.06660
 TITRANT H+ /M = 0.05019


| V/ML | E/MV | LG | NBAR | TH | TH-NBAR |
|------|------|-------|------|------|---------|
| 0.80 | -9.5 | -7.36 | 0.81 | 0.77 | -0.04 |
| 1.00 | 5.2 | -7.59 | 0.73 | 0.69 | -0.04 |
| 1.20 | 5.2 | -7.69 | 0.66 | 0.61 | -0.05 |
| 1.40 | 9.7 | -7.78 | 0.58 | 0.56 | -0.02 |
| 1.60 | 15.4 | -7.92 | 0.51 | 0.48 | -0.03 |
| 1.80 | 21.2 | -8.06 | 0.43 | 0.39 | -0.04 |
| 2.00 | 25.5 | -8.17 | 0.36 | 0.34 | -0.02 |
| 2.20 | 31.6 | -8.34 | 0.28 | 0.26 | -0.02 |
| 2.40 | 38.9 | -8.55 | 0.21 | 0.18 | -0.03 |
| 2.60 | 45.8 | -8.75 | 0.13 | 0.12 | -0.01 |

TITRATION 2

E0 /MV = 332.1
 INIT VOL /ML = 20.00
 H+ IN CELL/MMOL = 0.25100
 LIG IN CELL/MMOL = 0.09700
 M IN CELL /MMOL = 0.06600
 TITRANT B- /M = 0.05150

| | | | | | |
|------|-------|-------|------|------|-------|
| 4.26 | -51.2 | -3.53 | 0.32 | 0.31 | 0.00 |
| 4.30 | -54.9 | -3.50 | 0.33 | 0.33 | -0.01 |
| 4.40 | -58.9 | -3.49 | 0.35 | 0.33 | -0.01 |

Log K1(BaII)-DHP-15-aneN2O3

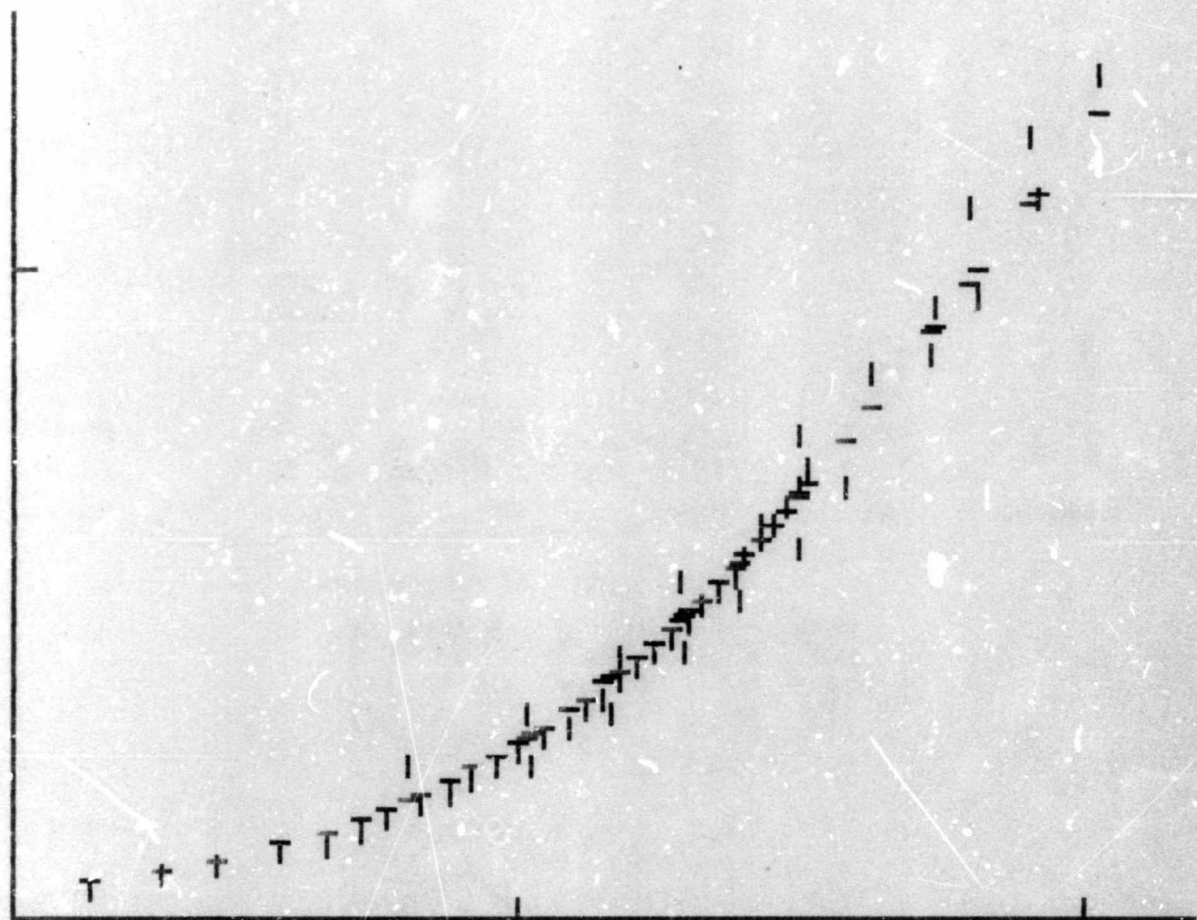

X - AXIS STARTS AT -4.9 UNITS = 1.0
 Y - AXIS STARTS AT 0 UNITS = 0.5

| = EXPERIMENTAL NBARS
 - = THEORETICAL NBARS

| | | | | | |
|------|-------|-------|------|------|-------|
| 4.26 | -51.2 | -3.53 | 0.32 | 0.31 | 0.00 |
| 4.30 | -54.9 | -3.50 | 0.33 | 0.33 | -0.01 |
| 4.40 | -58.9 | -3.49 | 0.35 | 0.33 | -0.01 |

AWT EWT 20 18.14 14 14.16.15

UNIVERSITY OF CALIFORNIA
 LAWRENCE BERKELEY LABORATORY
 CHEMICAL ENGINEERING DEPARTMENT
 BERKELEY, CALIFORNIA 94720

E0 /MV = 332.4
 INIT VOL /ML = 19.00
 LIG IN CELL/MMOL = 0.07298
 M IN CELL /MMOL = 0.06660
 TITRANT H+ /M = 0.05019

| V/ML | E/MV | LG | NBAR | TH | TH-NBAR |
|------|--------|-------|------|------|---------|
| 0.20 | -181.0 | -2.97 | 0.65 | 0.62 | -0.03 |
| 0.40 | -168.4 | -3.10 | 0.60 | 0.55 | -0.05 |
| 0.60 | -148.0 | -3.21 | 0.55 | 0.49 | -0.06 |
| 0.80 | -140.5 | -3.27 | 0.47 | 0.46 | -0.01 |
| 1.00 | -132.4 | -3.38 | 0.42 | 0.39 | -0.03 |
| 1.20 | -125.0 | -3.50 | 0.37 | 0.33 | -0.04 |
| 1.40 | -120.0 | -3.58 | 0.30 | 0.29 | -0.01 |
| 1.60 | -113.2 | -3.71 | 0.26 | 0.23 | -0.03 |
| 1.80 | -108.0 | -3.82 | 0.20 | 0.19 | -0.01 |
| 2.00 | -101.0 | -3.99 | 0.16 | 0.14 | -0.02 |
| 2.20 | -93.0 | -4.20 | 0.12 | 0.09 | -0.03 |

TITRATION 3

E0 /MV = 417.0
 INIT VOL /ML = 20.00
 H+ IN CELL/MMOL = 0.24900
 LIG IN CELL/MMOL = 0.09700
 M IN CELL /MMOL = 0.19980
 TITRANT B- /M = 0.05090

| V/ML | E/MV | LG | NBAR | TH | TH-NBAR |
|------|-------|-------|------|------|---------|
| 1.50 | 15.0 | -4.76 | 0.02 | 0.03 | 0.00 |
| 1.63 | 10.8 | -4.64 | 0.03 | 0.03 | 0.00 |
| 1.70 | 7.2 | -4.53 | 0.04 | 0.04 | 0.00 |
| 1.82 | 3.0 | -4.42 | 0.05 | 0.06 | 0.01 |
| 1.99 | 0.1 | -4.34 | 0.05 | 0.07 | 0.01 |
| 2.00 | -2.3 | -4.28 | 0.06 | 0.07 | 0.01 |
| 2.10 | -4.2 | -4.24 | 0.07 | 0.08 | 0.01 |
| 2.20 | -6.7 | -4.18 | 0.08 | 0.09 | 0.01 |
| 2.30 | -9.0 | -4.12 | 0.09 | 0.10 | 0.01 |
| 2.40 | -11.0 | -4.08 | 0.10 | 0.11 | 0.01 |
| 2.50 | -13.0 | -4.04 | 0.12 | 0.12 | 0.01 |
| 2.60 | -15.0 | -4.00 | 0.13 | 0.13 | 0.01 |
| 2.70 | -17.0 | -3.96 | 0.14 | 0.15 | 0.01 |
| 2.80 | -19.2 | -3.91 | 0.15 | 0.16 | 0.01 |
| 2.90 | -21.0 | -3.88 | 0.16 | 0.17 | 0.01 |
| 3.00 | -23.0 | -3.85 | 0.17 | 0.18 | 0.01 |
| 3.10 | -24.9 | -3.82 | 0.18 | 0.19 | 0.01 |
| 3.20 | -26.8 | -3.79 | 0.19 | 0.20 | 0.01 |
| 3.30 | -28.8 | -3.76 | 0.20 | 0.21 | 0.01 |
| 3.40 | -30.8 | -3.73 | 0.22 | 0.22 | 0.01 |
| 3.50 | -33.0 | -3.70 | 0.23 | 0.24 | 0.01 |
| 3.60 | -35.0 | -3.68 | 0.24 | 0.24 | 0.00 |
| 3.70 | -37.5 | -3.65 | 0.25 | 0.26 | 0.01 |
| 3.80 | -39.9 | -3.62 | 0.26 | 0.27 | 0.01 |
| 3.90 | -42.3 | -3.60 | 0.26 | 0.28 | 0.00 |
| 4.00 | -45.0 | -3.57 | 0.29 | 0.29 | 0.00 |
| 4.10 | -48.0 | -3.55 | 0.30 | 0.20 | 0.00 |

Log K1(BaII)-DHP-15-aneN203

ACID CONSTANTS:

PKA (1) = 8.30

PKA (2) = 7.62

MATRIX CONDITION NUMBER = 1.0000

CONVERGENCE IN *2* CYCLES

RMSD = .0207

SUM DELTA SQUARED = .02

PARAMETERS AND ERRORS

LOG K (1) = 3.188 +- .007

TITRATION 1

E0 /MV = 417.7
 INIT VOL /ML = 20.00
 H+ IN CELL/MMOL = 0.25100
 LIG IN CELL/MMOL = 0.09700
 M IN CELL /MMOL = 0.06660
 TITRANT B- /M = 0.05150

| V/ML | E/MV | LG | NBAR | TH | TH-NBAR |
|------|-------|-------|------|------|---------|
| 2.00 | -11.2 | -3.98 | 0.12 | 0.14 | 0.02 |
| 2.20 | -17.1 | -3.83 | 0.16 | 0.18 | 0.03 |
| 2.42 | -22.8 | -3.70 | 0.20 | 0.23 | 0.03 |
| 2.60 | -27.2 | -3.61 | 0.24 | 0.27 | 0.03 |
| 2.80 | -32.6 | -3.50 | 0.28 | 0.33 | 0.04 |
| 3.00 | -37.2 | -3.42 | 0.33 | 0.37 | 0.04 |
| 3.40 | -46.9 | -3.27 | 0.44 | 0.45 | 0.02 |
| 3.66 | -52.8 | -3.19 | 0.48 | 0.50 | 0.02 |
| 3.90 | -62.2 | -3.08 | 0.55 | 0.56 | 0.01 |


TITRATION 2

X - AXIS STARTS AT -5.1 UNITS = 1.0

Y - AXIS STARTS AT 0 UNITS = 0.5


! = EXPERIMENTAL NBARS

- = THEORETICAL NBARS

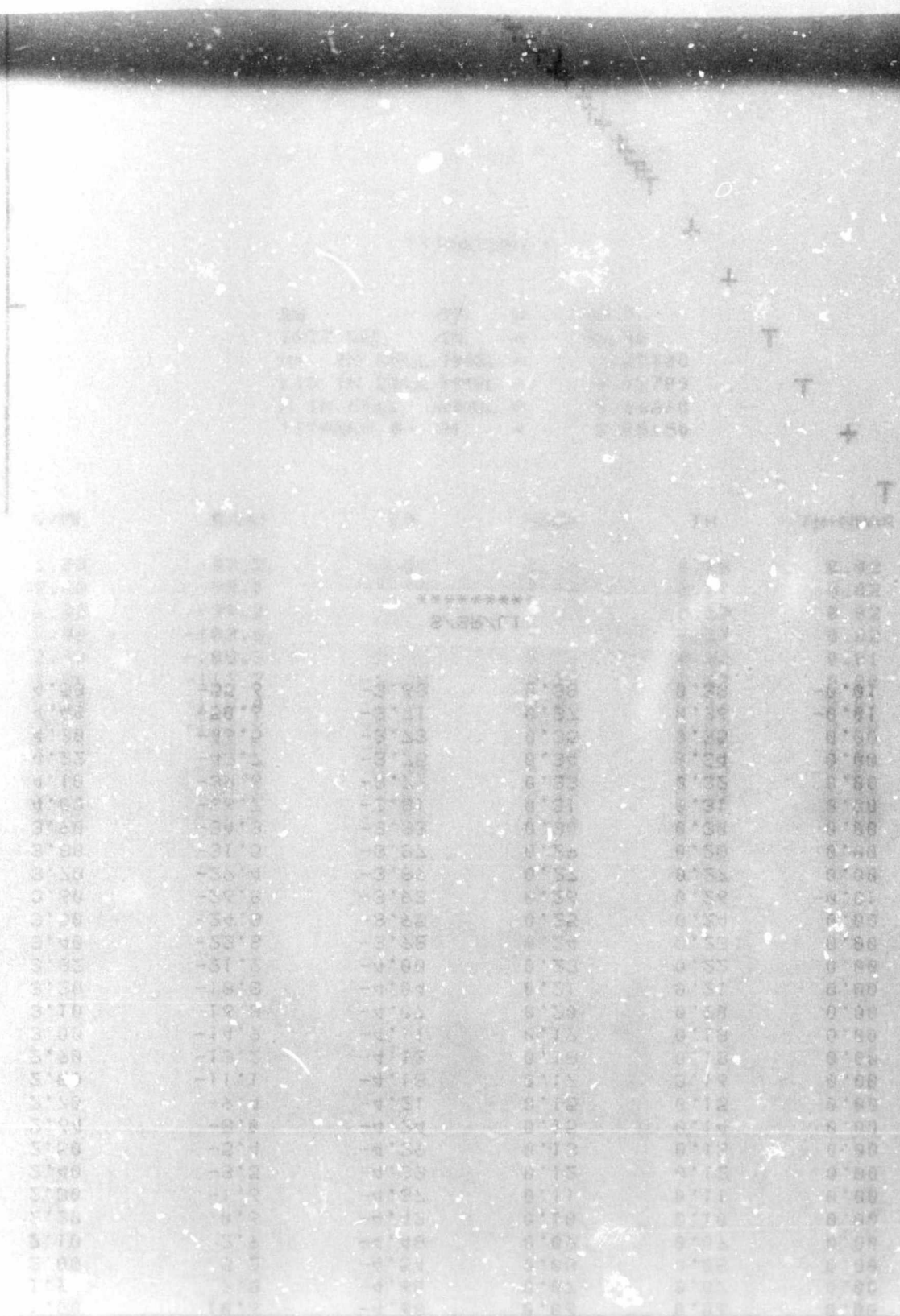

| | | | | | |
|------|-------|-------|------|------|-------|
| 1.80 | 10.6 | -4.68 | 0.06 | 0.06 | 0.00 |
| 1.90 | 7.8 | -4.60 | 0.07 | 0.07 | 0.00 |
| 2.00 | 5.2 | -4.54 | 0.08 | 0.08 | 0.00 |
| 2.10 | 2.9 | -4.48 | 0.09 | 0.09 | 0.00 |
| 2.20 | 0.6 | -4.42 | 0.10 | 0.10 | 0.00 |
| 2.30 | -1.6 | -4.37 | 0.11 | 0.11 | 0.00 |
| 2.40 | -3.5 | -4.33 | 0.12 | 0.12 | 0.00 |
| 2.50 | -5.4 | -4.29 | 0.13 | 0.13 | 0.00 |
| 2.64 | -8.0 | -4.24 | 0.15 | 0.14 | 0.00 |
| 2.70 | -9.4 | -4.21 | 0.15 | 0.15 | 0.00 |
| 2.80 | -11.1 | -4.18 | 0.17 | 0.16 | 0.00 |
| 2.90 | -13.2 | -4.13 | 0.18 | 0.18 | 0.00 |
| 3.00 | -14.9 | -4.11 | 0.19 | 0.18 | 0.00 |
| 3.10 | -16.8 | -4.07 | 0.20 | 0.20 | 0.00 |
| 3.20 | -18.8 | -4.04 | 0.21 | 0.21 | 0.00 |
| 3.32 | -21.2 | -4.00 | 0.23 | 0.22 | 0.00 |
| 3.40 | -22.8 | -3.98 | 0.24 | 0.23 | 0.00 |
| 3.50 | -24.8 | -3.95 | 0.25 | 0.24 | 0.00 |
| 3.60 | -26.8 | -3.92 | 0.26 | 0.26 | -0.01 |
| 3.70 | -29.4 | -3.89 | 0.27 | 0.27 | 0.00 |
| 3.80 | -31.5 | -3.87 | 0.29 | 0.28 | 0.00 |
| 3.90 | -34.3 | -3.83 | 0.30 | 0.30 | 0.00 |
| 4.00 | -36.7 | -3.81 | 0.31 | 0.31 | 0.00 |
| 4.10 | -39.6 | -3.79 | 0.33 | 0.32 | 0.00 |
| 4.22 | -43.7 | -3.75 | 0.34 | 0.34 | 0.00 |
| 4.30 | -46.5 | -3.73 | 0.35 | 0.35 | 0.00 |
| 4.40 | -50.6 | -3.71 | 0.37 | 0.36 | -0.01 |
| 4.50 | -55.6 | -3.68 | 0.38 | 0.38 | -0.01 |

Log K1(SrII)-DHP-15-aneN203

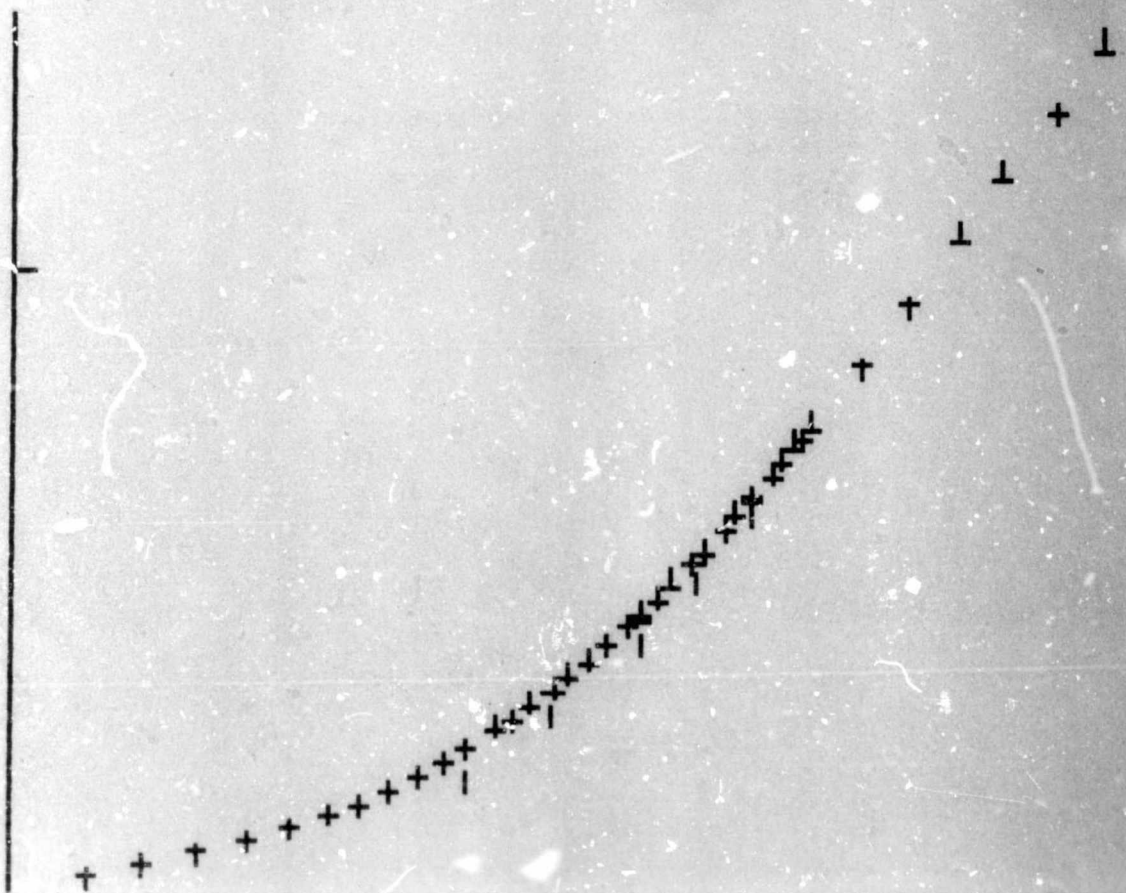

ACID CONSTANTS:

PKA (1) = 8.30

PKA (2) = 7.62

MATRIX CONDITION NUMBER = 1.0000

CONVERGENCE IN *2* CYCLES

RMSD = .0074

SUM DELTA SQUARED = .00

PARAMETERS AND ERRORS

LOG K (1) = 3.463 +- .003

TITRATION 1


E0 /MV = 330.7
 INIT VOL /ML = 20.00
 H+ IN CELL/MMOL = 0.25100
 LIG IN CELL/MMOL = 0.09700
 M IN CELL /MMOL = 0.06660
 TITRANT B- /M = 0.05150

| V/ML | E/MV | LG | NBAR | TH | TH-NBAR |
|------|--------|-------|------|------|---------|
| 1.80 | -87.2 | -4.29 | 0.11 | 0.13 | 0.02 |
| 2.00 | -93.0 | -4.14 | 0.16 | 0.17 | 0.02 |
| 2.22 | -99.2 | -3.99 | 0.21 | 0.23 | 0.02 |
| 2.40 | -103.6 | -3.89 | 0.26 | 0.27 | 0.02 |
| 2.60 | -108.3 | -3.78 | 0.31 | 0.32 | 0.01 |
| 2.80 | -112.7 | -3.70 | 0.37 | 0.37 | 0.00 |
| 3.02 | -118.4 | -3.59 | 0.42 | 0.43 | 0.01 |
| 3.20 | -122.9 | -3.51 | 0.47 | 0.47 | 0.00 |
| 3.42 | -128.3 | -3.42 | 0.53 | 0.52 | -0.01 |
| 3.60 | -133.8 | -3.34 | 0.58 | 0.57 | -0.01 |
| 3.80 | -140.6 | -3.25 | 0.62 | 0.62 | 0.00 |
| 4.02 | -148.4 | -3.16 | 0.68 | 0.67 | -0.01 |

TITRATION 2

E0 /MV = 417.1
 INIT VOL /ML = 21.00
 H+ IN CELL/MMOL = 0.24900
 LIG IN CELL/MMOL = 0.09700
 M IN CELL /MMOL = 0.19980
 TITRANT B- /M = 0.05000

| V/ML | E/MV | LG | NBAR | TH | TH-NBAR |
|------|------|-------|------|------|---------|
| 1.50 | 20.8 | -4.96 | 0.03 | 0.03 | 0.00 |
| 1.60 | 17.3 | -4.86 | 0.04 | 0.04 | 0.00 |
| 1.70 | 13.7 | -4.76 | 0.05 | 0.05 | 0.00 |

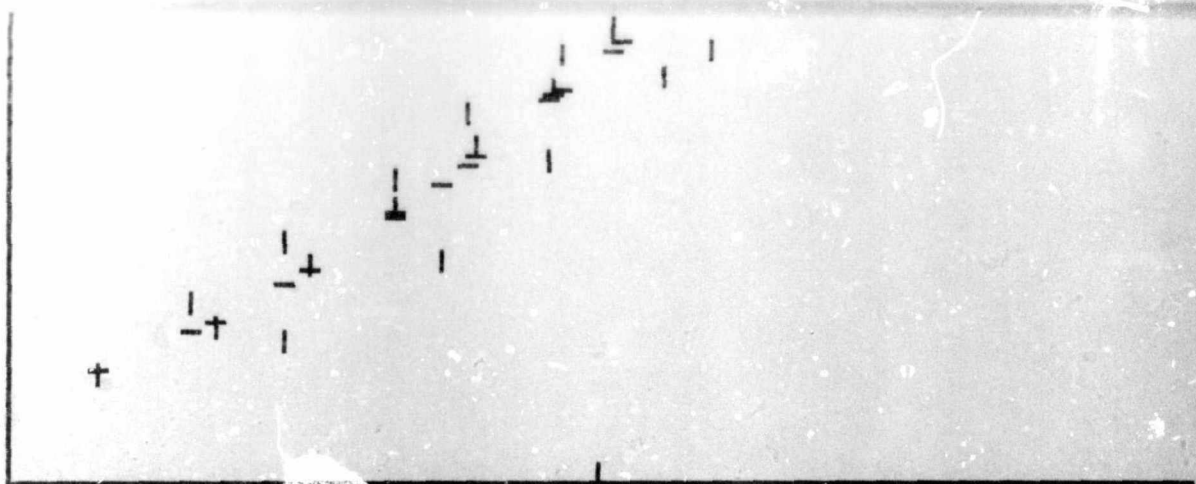

X - AXIS STARTS AT -1.9 UNITS = 1.0
 Y - AXIS STARTS AT 0 UNITS = 0.5

+ = EXPERIMENTAL NBARS
 - = THEORETICAL NBARS

| | | | | | |
|------|------|------|------|------|------|
| 0.00 | 0.00 | 0.00 | 0.00 | 0.00 | 0.00 |
| 0.05 | 0.05 | 0.05 | 0.05 | 0.05 | 0.05 |
| 0.10 | 0.10 | 0.10 | 0.10 | 0.10 | 0.10 |
| 0.15 | 0.15 | 0.15 | 0.15 | 0.15 | 0.15 |
| 0.20 | 0.20 | 0.20 | 0.20 | 0.20 | 0.20 |
| 0.25 | 0.25 | 0.25 | 0.25 | 0.25 | 0.25 |
| 0.30 | 0.30 | 0.30 | 0.30 | 0.30 | 0.30 |
| 0.35 | 0.35 | 0.35 | 0.35 | 0.35 | 0.35 |
| 0.40 | 0.40 | 0.40 | 0.40 | 0.40 | 0.40 |
| 0.45 | 0.45 | 0.45 | 0.45 | 0.45 | 0.45 |
| 0.50 | 0.50 | 0.50 | 0.50 | 0.50 | 0.50 |
| 0.55 | 0.55 | 0.55 | 0.55 | 0.55 | 0.55 |
| 0.60 | 0.60 | 0.60 | 0.60 | 0.60 | 0.60 |
| 0.65 | 0.65 | 0.65 | 0.65 | 0.65 | 0.65 |
| 0.70 | 0.70 | 0.70 | 0.70 | 0.70 | 0.70 |
| 0.75 | 0.75 | 0.75 | 0.75 | 0.75 | 0.75 |
| 0.80 | 0.80 | 0.80 | 0.80 | 0.80 | 0.80 |
| 0.85 | 0.85 | 0.85 | 0.85 | 0.85 | 0.85 |
| 0.90 | 0.90 | 0.90 | 0.90 | 0.90 | 0.90 |
| 0.95 | 0.95 | 0.95 | 0.95 | 0.95 | 0.95 |
| 1.00 | 1.00 | 1.00 | 1.00 | 1.00 | 1.00 |

111001 0.0000 0.0000
 111002 0.0000 0.0000
 111003 0.0000 0.0000
 111004 0.0000 0.0000
 111005 0.0000 0.0000

| | | | | | |
|------|------|------|------|------|------|
| 0.00 | 0.00 | 0.00 | 0.00 | 0.00 | 0.00 |
| 0.05 | 0.05 | 0.05 | 0.05 | 0.05 | 0.05 |
| 0.10 | 0.10 | 0.10 | 0.10 | 0.10 | 0.10 |
| 0.15 | 0.15 | 0.15 | 0.15 | 0.15 | 0.15 |
| 0.20 | 0.20 | 0.20 | 0.20 | 0.20 | 0.20 |
| 0.25 | 0.25 | 0.25 | 0.25 | 0.25 | 0.25 |
| 0.30 | 0.30 | 0.30 | 0.30 | 0.30 | 0.30 |
| 0.35 | 0.35 | 0.35 | 0.35 | 0.35 | 0.35 |
| 0.40 | 0.40 | 0.40 | 0.40 | 0.40 | 0.40 |
| 0.45 | 0.45 | 0.45 | 0.45 | 0.45 | 0.45 |
| 0.50 | 0.50 | 0.50 | 0.50 | 0.50 | 0.50 |
| 0.55 | 0.55 | 0.55 | 0.55 | 0.55 | 0.55 |
| 0.60 | 0.60 | 0.60 | 0.60 | 0.60 | 0.60 |
| 0.65 | 0.65 | 0.65 | 0.65 | 0.65 | 0.65 |
| 0.70 | 0.70 | 0.70 | 0.70 | 0.70 | 0.70 |
| 0.75 | 0.75 | 0.75 | 0.75 | 0.75 | 0.75 |
| 0.80 | 0.80 | 0.80 | 0.80 | 0.80 | 0.80 |
| 0.85 | 0.85 | 0.85 | 0.85 | 0.85 | 0.85 |
| 0.90 | 0.90 | 0.90 | 0.90 | 0.90 | 0.90 |
| 0.95 | 0.95 | 0.95 | 0.95 | 0.95 | 0.95 |
| 1.00 | 1.00 | 1.00 | 1.00 | 1.00 | 1.00 |

| V/ML | E/MV | LG | NBAR | TH | TH-NBAR |
|------|--------|-------|------|------|---------|
| 1.62 | -71.7 | -4.76 | 0.10 | 0.11 | 0.01 |
| 1.80 | -78.1 | -4.58 | 0.15 | 0.16 | 0.01 |
| 2.00 | -83.6 | -4.44 | 0.21 | 0.21 | 0.00 |
| 2.20 | -88.8 | -4.31 | 0.27 | 0.26 | -0.01 |
| 2.40 | -93.9 | -4.17 | 0.33 | 0.32 | -0.01 |
| 2.60 | -98.9 | -4.07 | 0.39 | 0.38 | -0.01 |
| 2.80 | -103.3 | -3.98 | 0.45 | 0.43 | -0.02 |
| 3.00 | -108.2 | -3.88 | 0.51 | 0.49 | -0.02 |
| 3.20 | -113.0 | -3.79 | 0.58 | 0.54 | -0.04 |
| 3.40 | -118.7 | -3.63 | 0.63 | 0.60 | -0.03 |
| 3.60 | -125.0 | -3.58 | 0.69 | 0.66 | -0.03 |
| 3.80 | -130.7 | -3.50 | 0.75 | 0.70 | -0.05 |
| 4.00 | -138.8 | -3.39 | 0.80 | 0.75 | -0.05 |

TITRATION 3

E0 /MV = 332.5
 INIT VOL /ML = 19.00
 LIG IN CELL/MMOL = 0.07298
 M IN CELL /MMOL = 0.06660
 TITRANT H+ /M = 0.05019

| V/ML | E/MV | LG | NBAR | TH | TH-NBAR |
|------|--------|-------|------|------|---------|
| 0.20 | -165.2 | -3.29 | 0.82 | 0.79 | -0.03 |
| 0.40 | -143.2 | -3.49 | 0.77 | 0.70 | -0.06 |
| 0.60 | -132.7 | -3.58 | 0.69 | 0.66 | -0.03 |
| 0.80 | -122.4 | -3.74 | 0.63 | 0.57 | -0.06 |
| 1.00 | -116.0 | -3.83 | 0.56 | 0.52 | -0.03 |
| 1.20 | -108.6 | -3.97 | 0.50 | 0.44 | -0.05 |
| 1.40 | -103.4 | -4.06 | 0.43 | 0.39 | -0.04 |
| 1.60 | -96.9 | -4.20 | 0.37 | 0.31 | -0.05 |
| 1.80 | -92.0 | -4.31 | 0.30 | 0.26 | -0.03 |
| 2.00 | -85.3 | -4.48 | 0.24 | 0.20 | -0.04 |
| 2.20 | -79.6 | -4.62 | 0.18 | 0.15 | -0.03 |

Log K1/CA/-DHP-15-aneN203

Log K1(CaII)-DHP-15-aneN203

ACID CONSTANTS:

PKA (1) = 8.30
 PKA (2) = 7.62

MATRIX CONDITION NUMBER = 1.0000
 CONVERGENCE IN *3* CYCLES
 RMSD = .0521
 SUM DELTA SQUARED = .09

PARAMETERS AND ERRORS

LOG K (1) = 3.864 +- .018

TITRATION 1

E0 /MV = 416.6
 INIT VOL /ML = 20.00
 H+ IN CELL/MMOL = 0.07530
 LIG IN CELL/MMOL = 0.02910
 M IN CELL /MMOL = 0.01650
 TITRANT B- /M = 0.05150

| V/ML | E/MV | LG | NBAR | TH | TH-NBAR |
|------|-------|-------|------|------|---------|
| 0.60 | -12.3 | -4.48 | 0.14 | 0.20 | 0.06 |
| 0.70 | -21.9 | -4.24 | 0.22 | 0.30 | 0.08 |
| 0.80 | -29.7 | -4.07 | 0.32 | 0.38 | 0.06 |
| 0.90 | -38.4 | -3.90 | 0.40 | 0.48 | 0.08 |
| 0.74 | -42.3 | -3.83 | 0.43 | 0.52 | 0.09 |
| 1.00 | -47.4 | -3.75 | 0.49 | 0.57 | 0.08 |
| 1.04 | -51.6 | -3.68 | 0.51 | 0.60 | 0.09 |
| 1.10 | -58.2 | -3.58 | 0.55 | 0.66 | 0.11 |
| 1.20 | -66.8 | -3.48 | 0.66 | 0.71 | 0.05 |

TITRATION 2

E0 /MV = 330.7
 INIT VOL /ML = 20.00
 H+ IN CELL/MMOL = 0.25100
 LIG IN CELL/MMOL = 0.09700
 M IN CELL /MMOL = 0.06650
 TITRANT B- /M = 0.05150

E0 /MV = 428.5
 INIT VOL /ML = 18.00
 H+ IN CELL/MMOL = 0.21795
 LIG IN CELL/MMOL = 0.07265
 M IN CELL /MMOL = 0.09990
 TITRANT B- /M = 0.05010

| V/ML | E/MV | LG | NBAR | TH | TH-NBAR |
|------|-------|--------|------|------|---------|
| 2.00 | 148.2 | -10.85 | 0.14 | 0.14 | 0.00 |
| 2.10 | 144.9 | -10.76 | 0.17 | 0.17 | 0.00 |
| 2.20 | 142.0 | -10.68 | 0.19 | 0.19 | 0.00 |
| 2.30 | 139.2 | -10.61 | 0.21 | 0.22 | 0.00 |
| 2.40 | 136.6 | -10.54 | 0.24 | 0.24 | 0.00 |
| 2.50 | 134.1 | -10.48 | 0.26 | 0.27 | 0.01 |
| 2.60 | 131.8 | -10.43 | 0.29 | 0.29 | 0.01 |
| 2.70 | 129.4 | -10.38 | 0.31 | 0.32 | 0.01 |
| 2.80 | 127.3 | -10.34 | 0.34 | 0.34 | 0.00 |
| 2.90 | 124.9 | -10.29 | 0.36 | 0.37 | 0.00 |
| 3.00 | 122.6 | -10.24 | 0.39 | 0.39 | 0.00 |
| 3.10 | 120.3 | -10.20 | 0.41 | 0.42 | 0.00 |
| 3.20 | 118.2 | -10.17 | 0.44 | 0.44 | 0.00 |
| 3.30 | 115.5 | -10.12 | 0.46 | 0.46 | 0.00 |
| 3.40 | 113.5 | -10.09 | 0.49 | 0.48 | -0.01 |
| 3.50 | 110.2 | -10.03 | 0.51 | 0.51 | 0.00 |
| 3.60 | 107.5 | -10.00 | 0.54 | 0.53 | -0.01 |
| 3.70 | 104.7 | -9.97 | 0.56 | 0.55 | -0.02 |
| 3.80 | 101.7 | -9.94 | 0.59 | 0.56 | -0.03 |
| 3.90 | 98.5 | -9.92 | 0.61 | 0.57 | -0.04 |

TITRATION 3

E0 /MV = 443.4
 INIT VOL /ML = 15.00
 H+ IN CELL/MMOL = 0.21795
 LIG IN CELL/MMOL = 0.07265
 M IN CELL /MMOL = 0.16650
 TITRANT B- /M = 0.05010

| V/ML | E/MV | LG | NBAR | TH | TH-NBAR |
|------|-------|--------|------|------|---------|
| 2.10 | 170.0 | -11.03 | 0.10 | 0.10 | 0.00 |
| 2.20 | 167.4 | -10.97 | 0.11 | 0.11 | -0.01 |
| 2.30 | 164.9 | -10.90 | 0.13 | 0.12 | -0.01 |
| 2.40 | 162.9 | -10.86 | 0.14 | 0.13 | -0.01 |
| 2.50 | 160.3 | -10.80 | 0.16 | 0.15 | -0.01 |
| 2.60 | 158.3 | -10.76 | 0.17 | 0.17 | -0.01 |
| 2.70 | 156.1 | -10.71 | 0.19 | 0.18 | -0.01 |
| 2.80 | 154.1 | -10.67 | 0.20 | 0.19 | -0.01 |
| 2.90 | 152.2 | -10.64 | 0.22 | 0.21 | -0.01 |
| 3.00 | 150.3 | -10.61 | 0.23 | 0.22 | -0.02 |
| 3.10 | 148.4 | -10.58 | 0.25 | 0.23 | -0.02 |
| 3.20 | 146.3 | -10.55 | 0.26 | 0.24 | -0.02 |
| 3.30 | 144.3 | -10.52 | 0.28 | 0.25 | -0.02 |
| 3.40 | 142.2 | -10.50 | 0.29 | 0.27 | -0.03 |
| 3.50 | 139.9 | -10.47 | 0.31 | 0.28 | -0.03 |

ACID CONSTANTS:

PKA (1) = 9.29
 PKA (2) = 8.50
 PKA (3) = 2.12

MATRIX CONDITION NUMBER = 1.0000
 CONVERGENCE IN *2* CYCLES
 RMSD = .0139
 SUM DELTA SQUARED = .01

PARAMETERS AND ERRORS


LOG K (1) = 10.053 +- .004

TITRATION 1

E0 /MV = 409.2
 INIT VOL /ML = 17.00
 H+ IN CELL/MMOL = 0.21795
 LIG IN CELL/MMOL = 0.07265
 M IN CELL /MMOL = 0.06660
 TITRANT B- /M = 0.05010

| U/ML | E/MV | LG | NBAR | TH | TH-NBAR |
|------|-------|--------|------|------|---------|
| 1.70 | 134.7 | -10.96 | 0.10 | 0.11 | 0.01 |
| 1.80 | 130.0 | -10.82 | 0.14 | 0.15 | 0.01 |
| 1.90 | 125.6 | -10.69 | 0.17 | 0.19 | 0.01 |
| 2.00 | 121.9 | -10.59 | 0.21 | 0.23 | 0.02 |
| 2.10 | 118.8 | -10.50 | 0.25 | 0.26 | 0.02 |
| 2.20 | 115.5 | -10.41 | 0.28 | 0.30 | 0.02 |
| 2.30 | 112.9 | -10.35 | 0.32 | 0.34 | 0.02 |
| 2.40 | 110.0 | -10.27 | 0.36 | 0.38 | 0.02 |
| 2.50 | 107.4 | -10.21 | 0.40 | 0.41 | 0.01 |
| 2.60 | 104.5 | -10.14 | 0.43 | 0.45 | 0.02 |
| 2.70 | 101.8 | -10.08 | 0.47 | 0.49 | 0.02 |
| 2.80 | 98.7 | -10.00 | 0.51 | 0.53 | 0.02 |
| 2.90 | 96.3 | -9.95 | 0.55 | 0.56 | 0.01 |
| 3.00 | 93.2 | -9.88 | 0.58 | 0.60 | 0.02 |
| 3.10 | 90.3 | -9.82 | 0.62 | 0.63 | 0.01 |
| 3.20 | 87.2 | -9.75 | 0.66 | 0.67 | 0.01 |
| 3.30 | 83.7 | -9.67 | 0.70 | 0.71 | 0.01 |
| 3.40 | 79.7 | -9.58 | 0.73 | 0.75 | 0.01 |
| 3.50 | 76.5 | -9.53 | 0.77 | 0.77 | 0.00 |
| 3.60 | 72.4 | -9.44 | 0.81 | 0.80 | -0.01 |
| 3.70 | 65.9 | -9.29 | 0.85 | 0.85 | 0.01 |
| 3.80 | 61.0 | -9.20 | 0.88 | 0.88 | -0.01 |
| 3.90 | 52.9 | -9.01 | 0.92 | 0.92 | 0.00 |
| 4.00 | 40.3 | -8.70 | 0.96 | 0.96 | 0.00 |

TITRATION 2

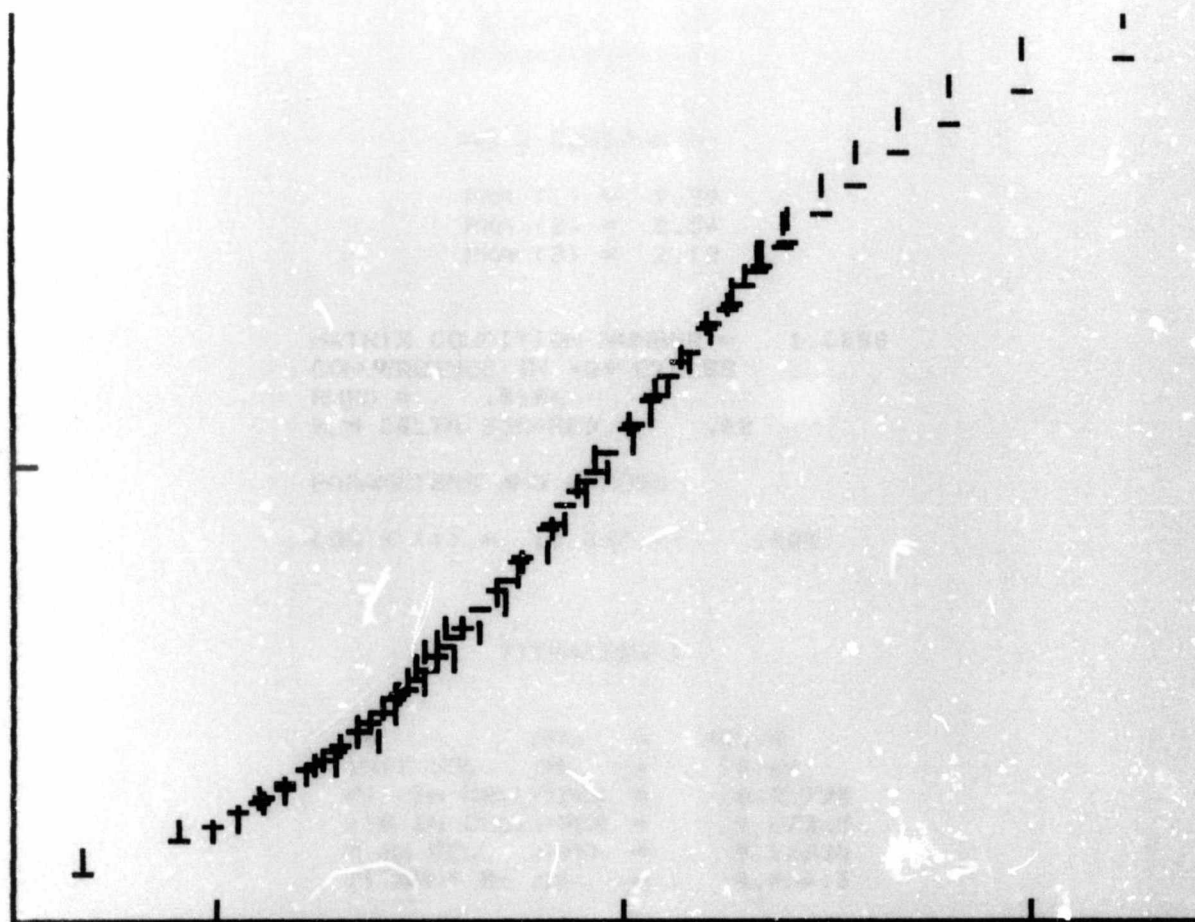

X - AXIS STARTS AT -11.5 UNITS = 1.0
 Y - AXIS STARTS AT 0 UNITS = 0.5

| = EXPERIMENTAL NBARS
 _ = THEORETICAL NBARS

11.50
 11.00
 10.50
 10.00
 9.50
 9.00
 8.50
 8.00
 7.50
 7.00
 6.50
 6.00
 5.50
 5.00
 4.50
 4.00
 3.50
 3.00
 2.50
 2.00
 1.50
 1.00
 0.50
 0.00

| V/ML | E/MV | LG | NBAR | TH | TH-NBAR |
|------|-------|--------|------|------|---------|
| 1.60 | 146.2 | -11.33 | 0.06 | 0.05 | -0.01 |
| 1.70 | 138.4 | -11.09 | 0.10 | 0.09 | -0.01 |
| 1.80 | 132.2 | -10.90 | 0.14 | 0.13 | -0.01 |
| 1.90 | 127.5 | -10.76 | 0.17 | 0.17 | 0.00 |
| 2.00 | 123.8 | -10.65 | 0.21 | 0.21 | 0.00 |
| 2.10 | 120.2 | -10.55 | 0.25 | 0.25 | 0.00 |
| 2.20 | 116.7 | -10.45 | 0.28 | 0.29 | 0.01 |
| 2.30 | 114.2 | -10.39 | 0.32 | 0.32 | 0.00 |
| 2.40 | 111.0 | -10.31 | 0.36 | 0.36 | 0.01 |
| 2.50 | 108.4 | -10.24 | 0.40 | 0.40 | 0.00 |
| 2.60 | 105.5 | -10.17 | 0.43 | 0.44 | 0.01 |
| 2.72 | 102.8 | -10.11 | 0.48 | 0.47 | -0.01 |
| 2.80 | 100.9 | -10.07 | 0.51 | 0.50 | -0.01 |
| 2.90 | 97.5 | -9.99 | 0.55 | 0.54 | 0.00 |
| 3.00 | 94.6 | -9.93 | 0.58 | 0.58 | 0.00 |
| 3.10 | 91.8 | -9.87 | 0.62 | 0.61 | -0.01 |
| 3.20 | 88.6 | -9.80 | 0.66 | 0.65 | -0.01 |
| 3.30 | 85.7 | -9.74 | 0.70 | 0.68 | -0.02 |
| 3.40 | 82.3 | -9.67 | 0.73 | 0.71 | -0.02 |
| 3.50 | 78.7 | -9.60 | 0.77 | 0.75 | -0.03 |
| 3.60 | 74.7 | -9.52 | 0.81 | 0.78 | -0.03 |
| 3.70 | 70.3 | -9.44 | 0.85 | 0.81 | -0.04 |
| 3.80 | 65.0 | -9.33 | 0.88 | 0.84 | -0.04 |
| 3.90 | 58.6 | -9.20 | 0.92 | 0.88 | -0.04 |
| 4.00 | 50.2 | -9.03 | 0.96 | 0.92 | -0.04 |
| 4.10 | 38.3 | -8.78 | 1.00 | 0.95 | -0.05 |

TITRATION 3

E_0 /MV = 409.4
 INIT VOL /ML = 20.00
 H^+ IN CELL/MMOL = 0.21795
 LIG IN CELL/MMOL = 0.07265
 M IN CELL /MMOL = 0.08325
 TITRANT B- /M = 0.05010

| V/ML | E/MV | LG | NBAR | TH | TH-NBAR |
|------|-------|--------|------|------|---------|
| 1.80 | 131.9 | -10.94 | 0.11 | 0.12 | 0.01 |
| 1.90 | 128.1 | -10.83 | 0.14 | 0.15 | 0.01 |
| 2.00 | 123.7 | -10.70 | 0.17 | 0.19 | 0.02 |
| 2.10 | 120.1 | -10.60 | 0.20 | 0.23 | 0.03 |
| 2.20 | 118.2 | -10.56 | 0.23 | 0.24 | 0.02 |
| 2.30 | 115.3 | -10.48 | 0.26 | 0.28 | 0.02 |
| 2.40 | 112.4 | -10.41 | 0.29 | 0.31 | 0.02 |
| 2.50 | 109.9 | -10.35 | 0.32 | 0.34 | 0.03 |
| 2.60 | 107.3 | -10.29 | 0.35 | 0.38 | 0.03 |
| 2.70 | 105.5 | -10.25 | 0.38 | 0.39 | 0.02 |
| 2.80 | 102.5 | -10.18 | 0.41 | 0.43 | 0.03 |
| 2.90 | 100.3 | -10.14 | 0.44 | 0.46 | 0.02 |
| 3.00 | 98.0 | -10.09 | 0.47 | 0.48 | 0.02 |
| 3.10 | 95.3 | -10.04 | 0.50 | 0.52 | 0.02 |
| 3.20 | 92.5 | -9.98 | 0.53 | 0.55 | 0.02 |
| 3.30 | 90.0 | -9.94 | 0.56 | 0.57 | 0.02 |
| 3.40 | 87.4 | -9.90 | 0.59 | 0.60 | 0.01 |
| 3.50 | 84.3 | -9.84 | 0.62 | 0.63 | 0.01 |
| 3.60 | 81.1 | -9.79 | 0.65 | 0.65 | 0.01 |
| 3.70 | 78.1 | -9.75 | 0.68 | 0.67 | 0.00 |
| 3.80 | 74.5 | -9.70 | 0.71 | 0.70 | -0.01 |
| 3.90 | 70.6 | -9.66 | 0.74 | 0.72 | -0.02 |
| 4.00 | 66.0 | -9.62 | 0.77 | 0.74 | -0.03 |

Log K1(PbII)-15-ane-N302

ACID CONSTANTS:

PKA (1) = 9.29
PKA (2) = 8.50
PKA (3) = 2.12

MATRIX CONDITION NUMBER = 1.0000
CONVERGENCE IN *3* CYCLES
RMSD = .0184
SUM DELTA SQUARED = .02

PARAMETERS AND ERRORS

LOG K (1) = 10.067 +- .005

TITRATION 1

E0 /MV = 443.4
INIT VOL /ML = 15.00
H+ IN CELL/MMOL = 0.21795
LIG IN CELL/MMOL = 0.07265
M IN CELL /MMOL = 0.16650
TITRANT B- /M = 0.05010


| V/ML | E/MV | LG | NBAR | TH | TH-NBAR |
|------|-------|--------|------|------|---------|
| 2.10 | 169.4 | -11.01 | 0.10 | 0.10 | 0.00 |
| 2.20 | 166.8 | -10.74 | 0.11 | 0.12 | 0.00 |
| 2.30 | 164.3 | -10.88 | 0.13 | 0.13 | 0.00 |
| 2.40 | 162.0 | -10.83 | 0.14 | 0.15 | 0.00 |
| 2.50 | 159.7 | -10.78 | 0.16 | 0.16 | 0.00 |
| 2.60 | 157.7 | -10.74 | 0.17 | 0.18 | 0.00 |
| 2.70 | 155.6 | -10.69 | 0.19 | 0.19 | 0.00 |
| 2.80 | 153.6 | -10.65 | 0.20 | 0.21 | 0.00 |
| 2.90 | 151.7 | -10.62 | 0.22 | 0.22 | 0.00 |
| 3.00 | 149.8 | -10.59 | 0.23 | 0.23 | 0.00 |
| 3.10 | 147.9 | -10.56 | 0.25 | 0.24 | -0.01 |
| 3.20 | 145.9 | -10.53 | 0.26 | 0.25 | -0.01 |
| 3.30 | 143.9 | -10.51 | 0.28 | 0.27 | -0.01 |
| 3.40 | 141.8 | -10.48 | 0.29 | 0.28 | -0.02 |
| 3.50 | 139.6 | -10.46 | 0.31 | 0.29 | -0.02 |
| 3.60 | 137.2 | -10.43 | 0.32 | 0.30 | -0.02 |

TITRATION 2

E0 /MV = 449.2
INIT VOL /ML = 17.00
H+ IN CELL/MMOL = 0.21795
LIG IN CELL/MMOL = 0.07265
M IN CELL /MMOL = 0.16650
TITRANT B- /M = 0.05010

| | | | | | |
|------|--------|------|------|------|-------|
| 7.00 | -144.6 | 9.35 | 0.61 | 0.56 | -0.05 |
| 7.10 | -146.7 | 9.39 | 0.58 | 0.53 | -0.05 |
| 7.20 | -148.7 | 9.42 | 0.54 | 0.50 | -0.04 |
| 7.30 | -150.7 | 9.46 | 0.51 | 0.47 | -0.04 |
| 7.40 | -152.8 | 9.49 | 0.48 | 0.45 | -0.03 |
| 7.50 | -154.4 | 9.52 | 0.44 | 0.43 | -0.01 |

PKA-15-ane-N302

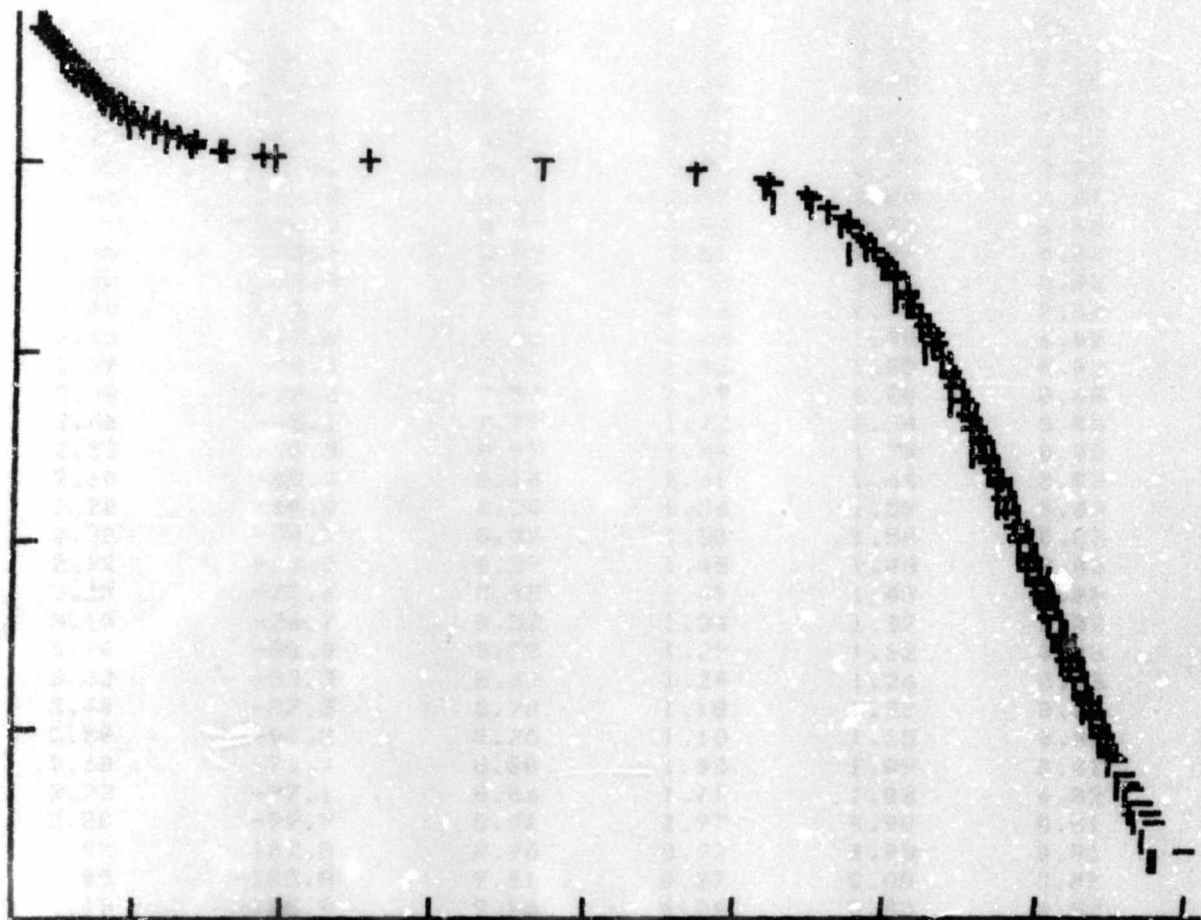

X - AXIS STARTS AT 2.2 UNITS = 1.0
 Y - AXIS STARTS AT 0 UNITS = 0.5

+ = EXPERIMENTAL NBARS
 - = THEORETICAL NBARS

| U/ML | E/MU | PH | NBAR | TH | TH-NBAR |
|------|--------|------|------|------|---------|
| 0.10 | 268.6 | 2.37 | 2.38 | 2.36 | -0.02 |
| 0.20 | 267.6 | 2.39 | 2.36 | 2.35 | -0.01 |
| 0.30 | 266.4 | 2.41 | 2.35 | 2.34 | -0.01 |
| 0.40 | 265.2 | 2.43 | 2.34 | 2.33 | -0.01 |
| 0.50 | 264.1 | 2.44 | 2.33 | 2.32 | -0.01 |
| 0.60 | 262.6 | 2.47 | 2.32 | 2.31 | -0.01 |
| 0.70 | 261.2 | 2.49 | 2.31 | 2.30 | -0.01 |
| 0.80 | 259.8 | 2.52 | 2.29 | 2.28 | -0.01 |
| 0.90 | 258.1 | 2.55 | 2.29 | 2.27 | -0.01 |
| 1.00 | 256.8 | 2.57 | 2.27 | 2.26 | -0.01 |
| 1.10 | 255.1 | 2.60 | 2.26 | 2.25 | -0.01 |
| 1.20 | 253.4 | 2.63 | 2.25 | 2.24 | -0.01 |
| 1.30 | 251.7 | 2.65 | 2.23 | 2.22 | -0.01 |
| 1.40 | 249.7 | 2.69 | 2.22 | 2.21 | -0.01 |
| 1.50 | 247.8 | 2.72 | 2.21 | 2.20 | -0.01 |
| 1.60 | 245.7 | 2.76 | 2.19 | 2.19 | -0.01 |
| 1.70 | 243.5 | 2.79 | 2.18 | 2.17 | -0.01 |
| 1.80 | 240.8 | 2.84 | 2.17 | 2.16 | -0.01 |
| 1.90 | 238.1 | 2.88 | 2.15 | 2.15 | -0.01 |
| 2.00 | 235.3 | 2.93 | 2.14 | 2.13 | -0.01 |
| 2.10 | 232.0 | 2.99 | 2.13 | 2.12 | -0.01 |
| 2.20 | 228.2 | 3.05 | 2.11 | 2.10 | -0.01 |
| 2.30 | 224.4 | 3.12 | 2.10 | 2.09 | 0.00 |
| 2.40 | 219.7 | 3.20 | 2.08 | 2.08 | 0.00 |
| 2.50 | 214.0 | 3.29 | 2.07 | 2.06 | 0.00 |
| 2.60 | 206.6 | 3.42 | 2.05 | 2.05 | 0.00 |
| 2.70 | 196.7 | 3.58 | 2.03 | 2.03 | 0.00 |
| 2.80 | 180.4 | 3.86 | 2.02 | 2.02 | 0.00 |
| 2.90 | 139.3 | 4.55 | 2.00 | 2.00 | 0.00 |
| 3.00 | 10.1 | 6.74 | 1.97 | 1.98 | 0.01 |
| 3.10 | -19.1 | 7.23 | 1.94 | 1.95 | 0.01 |
| 3.20 | -33.2 | 7.47 | 1.90 | 1.91 | 0.01 |
| 3.30 | -42.6 | 7.63 | 1.87 | 1.88 | 0.01 |
| 3.40 | -49.9 | 7.75 | 1.84 | 1.84 | 0.00 |
| 3.50 | -55.8 | 7.85 | 1.80 | 1.80 | 0.00 |
| 3.60 | -61.0 | 7.94 | 1.77 | 1.77 | 0.00 |
| 3.70 | -65.1 | 8.01 | 1.73 | 1.73 | 0.00 |
| 3.80 | -69.2 | 8.08 | 1.70 | 1.69 | 0.00 |
| 3.90 | -73.2 | 8.15 | 1.66 | 1.65 | -0.01 |
| 4.00 | -76.4 | 8.20 | 1.63 | 1.62 | -0.01 |
| 4.12 | -79.2 | 8.25 | 1.59 | 1.59 | 0.00 |
| 4.20 | -83.1 | 8.31 | 1.56 | 1.54 | -0.02 |
| 4.30 | -85.3 | 8.35 | 1.53 | 1.51 | -0.02 |
| 4.40 | -88.0 | 8.40 | 1.49 | 1.47 | -0.02 |
| 4.50 | -90.7 | 8.44 | 1.46 | 1.44 | -0.02 |
| 4.60 | -93.2 | 8.48 | 1.43 | 1.40 | -0.03 |
| 4.70 | -95.5 | 8.52 | 1.39 | 1.36 | -0.03 |
| 4.80 | -98.0 | 8.57 | 1.36 | 1.33 | -0.03 |
| 4.90 | -100.4 | 8.61 | 1.32 | 1.29 | -0.04 |
| 5.00 | -103.0 | 8.65 | 1.29 | 1.24 | -0.04 |
| 5.10 | -105.0 | 8.68 | 1.25 | 1.21 | -0.04 |
| 5.20 | -107.3 | 8.72 | 1.22 | 1.17 | -0.05 |
| 5.30 | -109.4 | 8.76 | 1.19 | 1.14 | -0.05 |
| 5.40 | -111.6 | 8.80 | 1.15 | 1.10 | -0.05 |
| 5.50 | -113.7 | 8.83 | 1.12 | 1.06 | -0.05 |
| 5.60 | -115.9 | 8.87 | 1.08 | 1.03 | -0.06 |
| 5.70 | -118.1 | 8.91 | 1.05 | 0.99 | -0.06 |
| 5.80 | -120.2 | 8.94 | 1.02 | 0.95 | -0.07 |
| 5.90 | -122.3 | 8.98 | 0.98 | 0.92 | -0.07 |
| 6.00 | -124.5 | 9.01 | 0.95 | 0.88 | -0.07 |
| 6.10 | -126.5 | 9.05 | 0.91 | 0.84 | -0.07 |
| 6.20 | -128.6 | 9.08 | 0.88 | 0.81 | -0.07 |
| 6.30 | -130.2 | 9.11 | 0.85 | 0.78 | -0.07 |
| 6.40 | -132.8 | 9.15 | 0.81 | 0.74 | -0.07 |
| 6.50 | -134.7 | 9.19 | 0.78 | 0.71 | -0.07 |
| 6.60 | -136.6 | 9.22 | 0.74 | 0.68 | -0.07 |
| 6.70 | -138.7 | 9.25 | 0.71 | 0.64 | -0.07 |
| 6.80 | -140.5 | 9.28 | 0.68 | 0.62 | -0.06 |
| 6.90 | -142.6 | 9.32 | 0.64 | 0.59 | -0.06 |

| V/ML | E/MV | PH | NBAR | TH | TH-NBAR |
|------|--------|------|------|------|---------|
| 0.10 | 279.3 | 2.50 | 2.27 | 2.29 | 0.02 |
| 0.20 | 277.6 | 2.53 | 2.26 | 2.28 | 0.02 |
| 0.30 | 275.9 | 2.56 | 2.24 | 2.27 | 0.02 |
| 0.40 | 274.0 | 2.59 | 2.23 | 2.25 | 0.02 |
| 0.50 | 272.0 | 2.63 | 2.21 | 2.24 | 0.02 |
| 0.60 | 269.8 | 2.66 | 2.20 | 2.22 | 0.02 |
| 0.80 | 265.0 | 2.74 | 2.17 | 2.19 | 0.02 |
| 0.90 | 262.3 | 2.79 | 2.16 | 2.18 | 0.02 |
| 1.00 | 259.3 | 2.84 | 2.14 | 2.16 | 0.02 |
| 1.10 | 256.1 | 2.89 | 2.12 | 2.14 | 0.02 |
| 1.22 | 252.3 | 2.96 | 2.10 | 2.13 | 0.03 |
| 1.30 | 248.1 | 3.03 | 2.09 | 2.11 | 0.02 |
| 1.40 | 243.1 | 3.11 | 2.08 | 2.09 | 0.02 |
| 1.50 | 237.2 | 3.21 | 2.06 | 2.07 | 0.02 |
| 1.60 | 227.8 | 3.37 | 2.05 | 2.05 | 0.01 |
| 1.70 | 214.9 | 3.59 | 2.03 | 2.03 | 0.00 |
| 1.80 | 193.6 | 3.95 | 2.01 | 2.01 | 0.00 |
| 1.90 | 68.8 | 5.72 | 1.98 | 2.00 | 0.02 |
| 2.00 | 0.9 | 7.21 | 1.93 | 1.95 | 0.02 |
| 2.10 | -17.0 | 7.51 | 1.88 | 1.90 | 0.03 |
| 2.20 | -28.9 | 7.71 | 1.82 | 1.85 | 0.03 |
| 2.30 | -38.2 | 7.87 | 1.77 | 1.80 | 0.03 |
| 2.40 | -45.1 | 7.99 | 1.72 | 1.74 | 0.03 |
| 2.50 | -50.3 | 8.07 | 1.66 | 1.70 | 0.03 |
| 2.60 | -55.6 | 8.16 | 1.61 | 1.64 | 0.03 |
| 2.70 | -60.4 | 8.24 | 1.56 | 1.59 | 0.03 |
| 2.80 | -64.7 | 8.32 | 1.50 | 1.54 | 0.03 |
| 2.90 | -68.8 | 8.39 | 1.45 | 1.48 | 0.03 |
| 3.00 | -72.6 | 8.45 | 1.40 | 1.43 | 0.03 |
| 3.10 | -76.7 | 8.52 | 1.34 | 1.37 | 0.02 |
| 3.20 | -80.0 | 8.58 | 1.29 | 1.32 | 0.03 |
| 3.30 | -83.3 | 8.63 | 1.24 | 1.26 | 0.03 |
| 3.40 | -87.3 | 8.70 | 1.18 | 1.20 | 0.01 |
| 3.50 | -90.3 | 8.75 | 1.13 | 1.15 | 0.02 |
| 3.60 | -93.4 | 8.80 | 1.08 | 1.09 | 0.02 |
| 3.72 | -97.1 | 8.86 | 1.01 | 1.03 | 0.02 |
| 3.80 | -99.7 | 8.91 | 0.97 | 0.98 | 0.01 |
| 3.90 | -102.8 | 8.96 | 0.92 | 0.93 | 0.01 |
| 4.00 | -105.8 | 9.01 | 0.87 | 0.88 | 0.01 |
| 4.10 | -108.8 | 9.06 | 0.81 | 0.83 | 0.02 |
| 4.20 | -111.4 | 9.11 | 0.76 | 0.79 | 0.03 |
| 4.30 | -115.6 | 9.18 | 0.71 | 0.72 | 0.01 |
| 4.40 | -118.5 | 9.23 | 0.65 | 0.67 | 0.02 |
| 4.50 | -121.4 | 9.27 | 0.60 | 0.63 | 0.02 |
| 4.62 | -124.7 | 9.33 | 0.54 | 0.58 | 0.04 |
| 4.70 | -127.7 | 9.38 | 0.50 | 0.53 | 0.04 |
| 4.80 | -130.9 | 9.44 | 0.45 | 0.49 | 0.04 |
| 4.90 | -134.0 | 9.49 | 0.39 | 0.45 | 0.06 |
| 5.02 | -137.8 | 9.55 | 0.33 | 0.40 | 0.07 |
| 5.10 | -141.0 | 9.61 | 0.29 | 0.37 | 0.08 |
| 5.20 | -144.3 | 9.66 | 0.24 | 0.33 | 0.09 |
| 5.30 | -147.8 | 9.72 | 0.19 | 0.30 | 0.11 |
| 5.40 | -151.6 | 9.79 | 0.14 | 0.26 | 0.12 |

TITRATION 3

E0 /MV = 488.7
 INIT VOL /ML = 20.00
 H+ IN CELL/MMOL = 0.43890
 LIC IN CELL/MMOL = 0.14420
 TITRANT B- /M = 0.05010

PKA-15-ane-N302

MATRIX CONDITION NUMBER = 3.0159
 CONVERGENCE IN *3* CYCLES
 RMSD = .0421
 SUM DELTA SQUARED = .26

PARAMETERS AND ERRORS

PKA (1) = 9.291 +- .013
 PKA (2) = 8.496 +- .013
 PKA (3) = 2.116 +- .017

TITRATION 1

E0 /MV = 409.4
 INIT VOL /ML = 20.00
 H+ IN CELL/MMOL = 0.21795
 LIG IN CELL/MMOL = 0.07265
 TITRANT B- /M = 0.05010

| V/ML | F/MV | PH | NBAR | TH | TH-NBAR |
|------|--------|-------|------|------|---------|
| 0.20 | 252.7 | 2.65 | 2.24 | 2.23 | -0.01 |
| 0.40 | 248.3 | 2.72 | 2.19 | 2.20 | 0.01 |
| 0.60 | 242.0 | 2.83 | 2.17 | 2.16 | 0.00 |
| 0.80 | 235.2 | 2.94 | 2.12 | 2.13 | 0.01 |
| 1.00 | 229.9 | 3.07 | 2.10 | 2.09 | -0.01 |
| 1.20 | 206.1 | 3.40 | 2.06 | 2.05 | -0.01 |
| 1.40 | 176.5 | 3.94 | 2.00 | 2.01 | 0.01 |
| 1.60 | -20.5 | 7.27 | 1.90 | 1.94 | 0.05 |
| 1.80 | -50.6 | 7.78 | 1.76 | 1.83 | 0.07 |
| 2.00 | -68.2 | 8.07 | 1.62 | 1.70 | 0.08 |
| 2.20 | -80.1 | 8.27 | 1.48 | 1.57 | 0.08 |
| 2.40 | -90.6 | 8.45 | 1.35 | 1.43 | 0.08 |
| 2.60 | -99.3 | 8.60 | 1.21 | 1.29 | 0.08 |
| 2.80 | -108.9 | 8.76 | 1.07 | 1.13 | 0.06 |
| 3.00 | -116.6 | 8.89 | 0.94 | 1.00 | 0.07 |
| 3.20 | -126.6 | 9.06 | 0.80 | 0.83 | 0.03 |
| 3.40 | -134.5 | 9.19 | 0.66 | 0.70 | 0.04 |
| 3.60 | -142.5 | 9.33 | 0.53 | 0.58 | 0.05 |
| 3.80 | -152.8 | 9.50 | 0.40 | 0.44 | 0.04 |
| 4.00 | -161.3 | 9.65 | 0.27 | 0.34 | 0.08 |
| 4.20 | -171.2 | 9.81 | 0.14 | 0.25 | 0.11 |
| 4.40 | -182.0 | 10.00 | 0.02 | 0.17 | 0.15 |

TITRATION 2

E0 /MV = 427.3
 INIT VOL /ML = 20.00
 H+ IN CELL/MMOL = 0.28107
 LIG IN CELL/MMOL = 0.09369
 TITRANT B- /M = 0.05010

Author Bhavan Rekha

Name of thesis Ligand Design Using The Neutral Oxygen Donor To Control Metal Ion Size Selectivity. 1989

PUBLISHER:

University of the Witwatersrand, Johannesburg

©2013

LEGAL NOTICES:

Copyright Notice: All materials on the University of the Witwatersrand, Johannesburg Library website are protected by South African copyright law and may not be distributed, transmitted, displayed, or otherwise published in any format, without the prior written permission of the copyright owner.

Disclaimer and Terms of Use: Provided that you maintain all copyright and other notices contained therein, you may download material (one machine readable copy and one print copy per page) for your personal and/or educational non-commercial use only.

The University of the Witwatersrand, Johannesburg, is not responsible for any errors or omissions and excludes any and all liability for any errors in or omissions from the information on the Library website.